



FEATURES:

Guardian Bassinets feature KYDEX® on all horizontal surfaces for durability, including top, base panel and chart board. KYDEX® is 4x thicker than standard thermofilm giving it superior impact strength and durability. KYDEX is certified GREENGUARD® Gold for Children and Schools. It also withstands a wide range of cleaners and disinfectants.

The PETG plastic tub has 4x greater impact resistance than acrylic tubs and has hooks for Trendelenburg positioning during medical procedures.

Stainless steel push handles allow greater view of infant. The push handle design eliminates the hard to get to areas to make all surface areas accessible and easy to clean.

Full extension drawer slides with anti-tip mechanism.

Drawer liners are removable for easy cleaning.

The base panel was designed to reduce the chances of caregivers hitting their shins while the bassinet is in motion and also to reduce the chance of siblings using the base as a step to get to the newborn.

Two 5" locking and two non-locking swivel poly casters.

Ships in 100% recyclable packaging. Manufactured in the USA.

SPECIFICATIONS

Handle height:	36in (92 cm)
Surface height:	30in (76 cm)
Cabinet width:	34in (86 cm)
Cabinet depth:	21.75in (53 cm)
Weight:	112lbs (51 kg)
Hardware options:	Black inset pulls or stainless steel bow pulls



Optional chart holder.




Optional pull-out chart board.



Shin saving base design.



All Guardian Bassinets come standard with removable plastic drawer liners for ease of cleaning. This may aid in infection prevention. The 6in deep drawer liners are made of high-impact ABS plastic and are field replaceable.

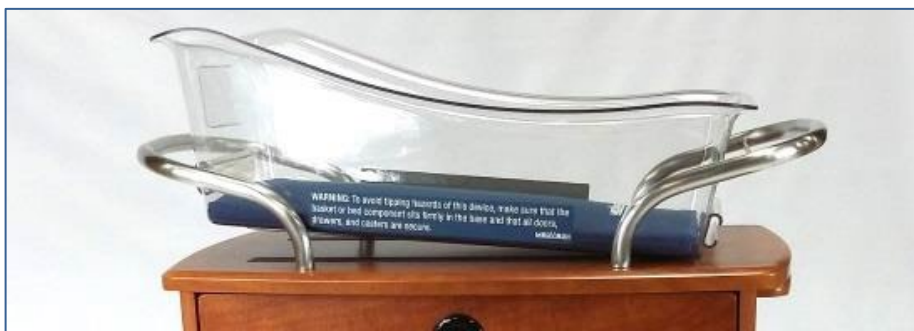

The Bassinet's tub and handle design allow for quick Trendelenburg or Reverse-Trendelenburg positioning during medical procedures.